



SINCE 1978

OPTICS TECHNOLOGY

Manufacturer & Exporters of : LAB TESTING INSTRUMENTS




An ISO 9001:2015 CERTIFIED IEC 61010-1:2010
An ISO 10002:2018 CERTIFIED BS EN 12469:2000
An ISO 13485:2016 CERTIFIED ISO 29463-5:2022



GLASS DISTILLATION (DOUBLE) CE

WITH AUTOMATIC WATER LEVEL CUT OFF

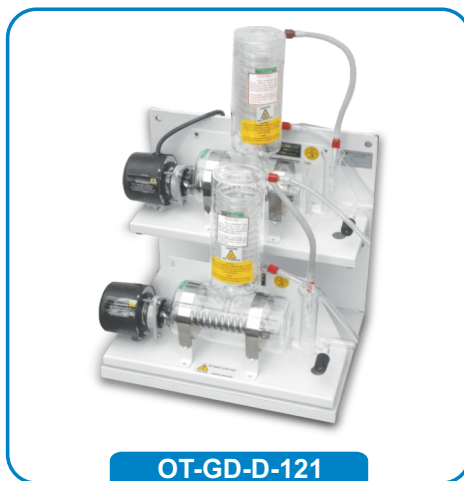

OT-GD-D-121

DOUBLE DISTILLATION

The distillation apparatus consists of boiler with heating elements embedded in glass and fused in spiral type coil internally of the bottom and tapered round glass, suitable to work on 220 Volts 50 Cycles AC complete with metallic stand rings and clamps to hold boiler condenser electric board with switch and connecting cords. It is capable of producing pure distilled water. **Also the whole unit can be hanged on wall.**

CAPACITY : 4 Ltr./hr. - Borosilicate Glass

2/3, DSIDC, (Behind Fun Cinema) NSP, Pitampura, Delhi-110034 (INDIA)

☎ 9811238150, 9811238151, 9811236960

E-mail: info@opticstechnology.in

opticstechnology@gmail.com

Website : www.opticstechnology.in