



SINCE 1978

# OPTICS TECHNOLOGY

Manufacturer & Exporters of : LAB TESTING INSTRUMENTS



An ISO 9001:2015 CERTIFIED IEC 61010-1:2010  
An ISO 10002:2018 CERTIFIED BS EN 12469:2000  
An ISO 13485:2016 CERTIFIED ISO 29463-5:2022




## WATER BATH INCUBATOR SHAKER CE

Heating Only (Metabolic Shaking Incubator)



OT-WBIS-60 MS

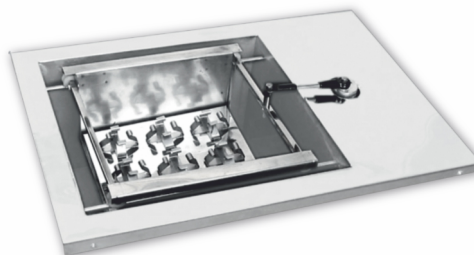

(SHAKING TRAY)

- Temperature range from room temperature to 100°C.
- Temp. is Controlled by Thermostat OR Microprocessor based PID Controller
- Double walled inside made of stainless steel 304 Quality and outside Mild Steel sheet painted in epoxy powder coating/S.S. 304.
- Oscillating tray is riding on ball rollers.
- Shaking speed range 40 to 140 RPM.
- To work on 220/230 volts A.C. Supply.

### SHAKING TRAY SIZE INSIDE

TRAY MADE OF S.S. 304 Quality

Model No.	Size
OT-WBIS-60-MS/SS-4	4 Flask of 250 ml.
OT-WBIS-60-MS/SS-6	6 Flask of 250 ml.
OT-WBIS-60-MS/SS-8	8 Flask of 250 ml.
OT-WBIS-60-MS/SS-12	12 Flask of 250 ml.

### Optional:

1. Microprocessor based PID Controller
2. Digital RPM indicator for speed
3. Water level cut off system.

MOC : MS-Mild steel duly powder coated SS-stainless steel kindly select the one required in model number.

2/3, DSIDC, (Behind Fun Cinema) NSP, Pitampura, Delhi-110034 (INDIA)

☎ 9811238150, 9811238151, 9811236960

E-mail: [info@opticstechnology.in](mailto:info@opticstechnology.in),

[opticstechnology@gmail.com](mailto:opticstechnology@gmail.com),

Website : [www.opticstechnology.in](http://www.opticstechnology.in)